

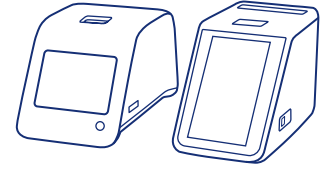




QuikRead go Strep A

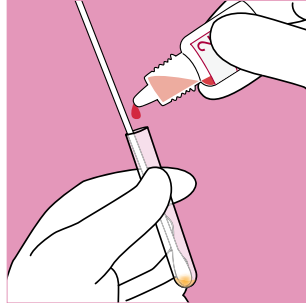



For QuikRead go Instrument and QuikRead go Plus Instrument

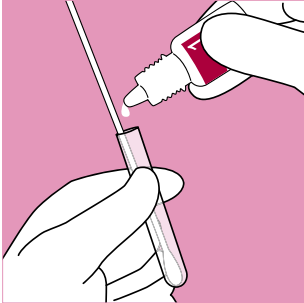



1 Take the throat sample with the swab intended for the kit. Place the swab into the extraction tube and add two drops of colourless Extraction reagent 1. Add two drops of the Extraction reagent 2. The solution turns yellow-orange.

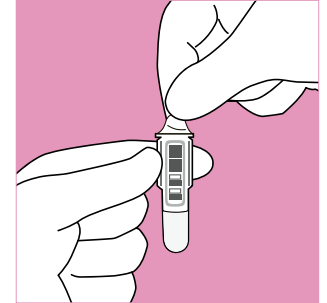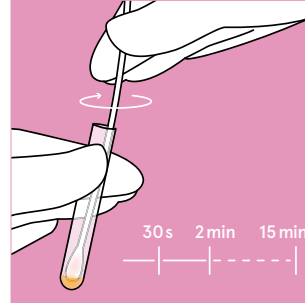


2 Swirl the swab in the solution for 30 s. Let it stand in the solution for an additional 1.5 minutes but no longer than a total of 15 minutes. Remove the foil cover from the cuvette. Do not touch the clear flat surfaces in the lower part of the cuvette.

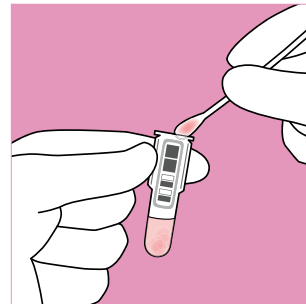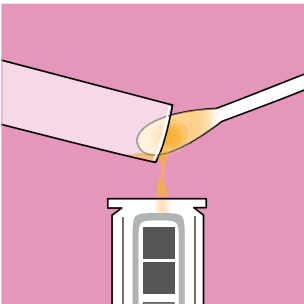

3 Tilt the extraction tube and slide the swab against its inner wall to pull the solution with the swab into the cuvette. If a large volume of liquid remains in the tube, pour it into the cuvette. Stir vigorously with the swab, the solution turns red. Press the swab against the inner wall when lifting it out.

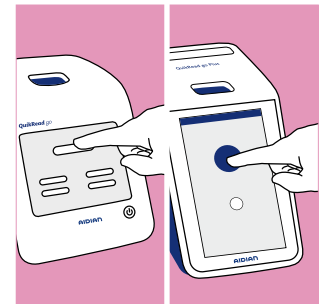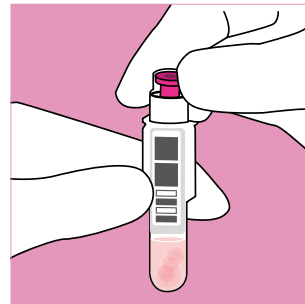


4 Close the cuvette tightly with the Strep A reagent cap. Do not press down the pink inner part. Perform the test within 4 hours. Choose Measure / Measure sample on the display of the QuikRead go instrument.

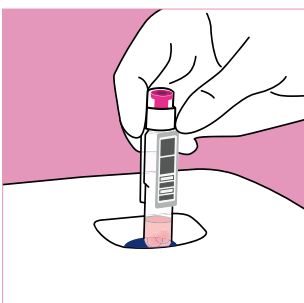

Please consult the instructions for use before performing the test.

5 Insert the cuvette into the measurement well of the instrument. The barcode should be facing you. Measurement time is less than 4 minutes. The result is displayed on the screen, and the cuvette rises automatically when the measurement is completed.

QuikRead go Strep A

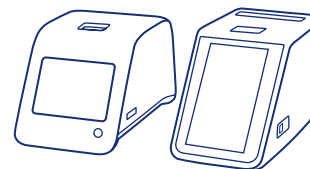

Kit

Unopened kit at 2–25°C: until the expiry date of the kit.



Reagent cap tube, extraction reagents and controls

Unopened kit components at 2–25°C: until the expiry date of the kit.
After the first opening at 2–25°C: 12 months.



Cuvettes

Unopened foil pouch at 2–25°C: until the expiry date of the kit.
Unopened cuvettes after opening the foil pouch:
at 2–8°C: 6 months / at 18–25°C: 3 months.




Opened cuvette


Let the buffer warm to 18–25°C before testing.
Begin the test within 2 hours from opening the cuvette.

Sampling instructions


When performing the test, swirling the swab during extraction and in the cuvette is important. Remember to press the swab against the inner wall when removing from cuvette.

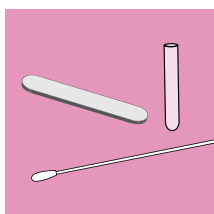


1 Collect everything you need for sampling close to the patient: HydraFlock flocked swab, tongue spatula and QuikRead go Strep A Extraction tube. Wear protected clothing including disposable gloves. Use only HydraFlock flocked swabs intended for the QuikRead go Strep A kit. Use of other swabs may yield false results, either positive or negative.

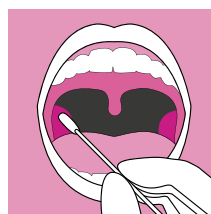

3 Collect the mucus/exudate with the swab. To ensure a good sample rotate the swab when sampling. Strep A infects tonsilles, remember to swab on both sides. Sample collected elsewhere may yield a negative result. A large amount of saliva or microbes from other parts of the mouth may disturb the test.

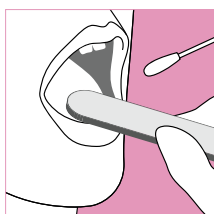

2 Remove the sterile HydraFlock flocked swab from its pouch. Ask the patient to open her/his mouth. Use the spatula to press down the tongue and start sampling.

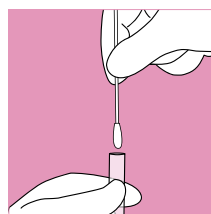

4 Take the swab out of the mouth and put it into the extraction tube. Be careful not to touch other surfaces of mouth to avoid contamination with normal flora.

Ensure correct extraction. Observe colour changes. The colour indicates a correct pH both during extraction and neutralization. Follow the given timelines.

Contact and orders

Aidian Oy, Finland

aidian@aidian.eu / +358 10 309 3000

orders@aidian.eu

quikread.com

Product	Cat. no.
QuikRead go Strep A	135883
QuikRead go Instrument	133893
QuikRead go Plus Instrument	155378